Provision of Guarantee

- 1 Any defects found within the warranty period under normal operation in compliance with the instruction manual shall be fixed and repaired free of charge.
 - 2 When requesting for free repair service, bring the product to the place of purchase, along with the warranty certificate.
 - 3 Warranty period is one year from the original purchase date.
 - 4 Even within the warranty period, repair service shall be performed at the customer's own cost in the following cases.
 - (A) Damage or malfunction from improper use
 - (B) Failure or damage due to careless handling such as falling after purchasing
 - (C) Damage or malfunction from remodeling or unauthorized repair
 - (D) Damage or malfunction from fire, earthquake, flood, lightning and other natural disasters or force majeure events
 - (E) If there is no presentation of this warranty
 - (F) If there is no serial NO, date of purchase, retailer, or customer name on this warranty, or If you rewrite the phrase.
 - (G) Accessories except the main body.
 - 5 These terms and conditions are valid only in Japan.
 - 6 This certificate shall not be reissued. Please keep in a secure place.
- ※ This certificate guarantees free repair service under the specified terms and conditions, and does not limit the customer's legal rights against the issuer or other businesses.
- ※ Any inquiries regarding the product or repair service after the warranty period shall be addressed to the place of purchase or sales agency.

Warranty Certificate	
Name of Product	N-LED FLEX
Serial No.	
Date of Purchase	
Warranty Period	One year from purchase (only applicable to main equipment)
Place of purchase	Store name · address · phone number Seal or Signature
Customer	Name
	Address 〒 Phone No.

Thank you for purchasing our product. This certificate guarantees the quality of the product for safe use. If any manufacturing defects are found within the warranty period under normal operation, please submit this warranty certificate to receive repair services free of charge.

- ※ If there is no store name and seal, date of purchase, or receipt, it will be invalid.
 ※ Date of purchase shall be entered by the store.

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ENGLISH

Instruction manual guidelines 10
 For safe and proper use 10
 Precautions before use 11
 Precautions for use 12
 Precautions for repair and inspection 13
 How to use 14
 Provision of Guarantee 31

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Aderans

Instruction manual guidelines

- This instruction manual is necessary for safe use of the product.
Please forward this Instruction Manual when you give this product to others.
- Read this instruction manual carefully before use.
- Keep in an easily accessible place for further reference.
- When abnormality occurs, please stop using it immediately and ask for repair.

For safe and proper use

○ Read the instructions carefully before use.

Pay extra attention to notices with warning signs.
Consequences of improper handling are classified as follows:

- | | | |
|---|----------------|--|
|  | Danger | Risk of death or serious injury |
|  | Warning | Possibility of death or serious injury |
|  | Caution | Possibility of personal or material damage |

 「Danger·Warning·Caution」requirements
※sign means「Risk of electric shock」

 「Restriction」which should be prohibited
※sign means「Do not dismantle」

 Which is a「Compulsory order」
※sign means「Pull the plug out」

Danger

-  Restriction
- Do not use in combination with the following medical devices:
 - Implantation devices such as pacemaker
 - Life support devices such as heart-lung machine
 - Installation devices such as electrocardiograph

Warning

●Do not use under the following conditions. This product is for beauty treatment only, and should not be used for any other purpose.

- Fever, sudden illness, infectious disease, heart disease, abnormal blood pressure
- Those pregnant or with malignant tumor
- Allergies (atopic dermatitis, etc.)
- Body parts that have undergone plastic surgery or with metal, plastic or silicon inside
- Those with hemophilia or hemorrhagic diseases
- Facial melanosis, freckles, warts
- Infection due to excessive sun exposure
- Photosensitive patients
- Infection from steroid hormones or cosmetics
- Children aged 12 and under

Restriction

Caution

- Do not use on wounds.
- Do not use under poor physical condition.
- Do not use after drinking alcohol.
- Consult a physician if you have an itchy scalp.
- Consult a physician before use if you have experienced any of the following:Restriction

- Allergic reaction to synthetic fiber (infection, etc.)
- Allergic reaction to cosmetics (infection, etc.)

Restriction

Precautions before use

Danger

- Do not use in combination with other beauty care devices.

Restriction

Warning

- Use only on designated parts. Do not use for purposes other than the original purpose.
- Those who are paralyzed or unable to communicate should not use this product.
- Do not look into the light when it is on.
- Do not use in wet or humid places.
- Do not place heavy objects on the AC adapter.
- Do not use power/voltage other than described in the specification.

Caution

- Fully charge before first use.
If the battery is not fully charged, the power may go out during use.
- Do not use in areas with high humidity and dust, as it may cause malfunction.
- If unused for a long time, check the normal and safe operation of the product before use.
- Do not use multiple wire connectors.
- Do not sprinkle or spill water on the product, as water may cause malfunction.
- Hold the plug head when plugging and unplugging the AC adapter. Pulling by the cord may cause snapping or disconnection.
- Do not unplug the AC adapter with wet hands, as it may cause electric shock.
- Remove accessories (earrings, hairpin, necklace, etc.)
- Do not use parts of other products.

Precautions for use

- Use properly according to the instructions.
Please conform to the precautions in the order of handling.

Warning

- In case abnormal symptoms or skin problems occur, stop using immediately and consult a physician.
- If excessive stimulation is felt, stop using immediately and consult a doctor.
- In the event of an earthquake or thunderstorm, turn the power off and unplug the AC adapter.
- Do not look into the light source while the LED light is on.

Restriction

⚠ Caution

Restriction

- If any defects are found, stop using immediately and inquire the store.
- Do not place the product near other electrical appliances (e.g. television, radio, PC, mobile phone, stereo, beauty care devices)
- Be careful not to drop the product during use.
- Do not use AC adapters other than the one included.
- Product may get hot during use, but it is not a safety issue.
- Product may fall down if the balance is tilted too much to one side.
- Do not use while connected to the AC adapter, as it may cause malfunction.

Precautions for repair and inspection

⚠ Warning

Restriction

- Do not dismantle, repair or remodel the product. In case of malfunctioning, stop using and inquire the store.
- Do not drop from a high place or give strong impact. In case of dropping by mistake, inspect the LED light and check for cracks before use.

⚠ Caution

Restriction


- Do not store in high temperature, humidity and dust level.
- Keep in a dry place out of reach of children.
- Do not store in a place with strong electromagnetic force.
- Do not damage, bend, pull or tie the AC adapter. Placing heavy objects on the AC adapter may damage the cord, and may cause fire or electric shock.
- Wash the outer cover according to the following instructions:
 - Hand-wash with neutral detergent. (Strong rubbing, squeezing or using the washing machine may damage the cover.)
 - Dry in a shady place with good ventilation.
 - The outer cover is not for permanent use. When damaged by frequent washing, purchase a new cover.

Regular inspection shall be performed.

- If unused for more than 6 months, check the normal and safe operation of the product before use. Using a malfunctioning product may lead to accidents.

How to use

Name of each part

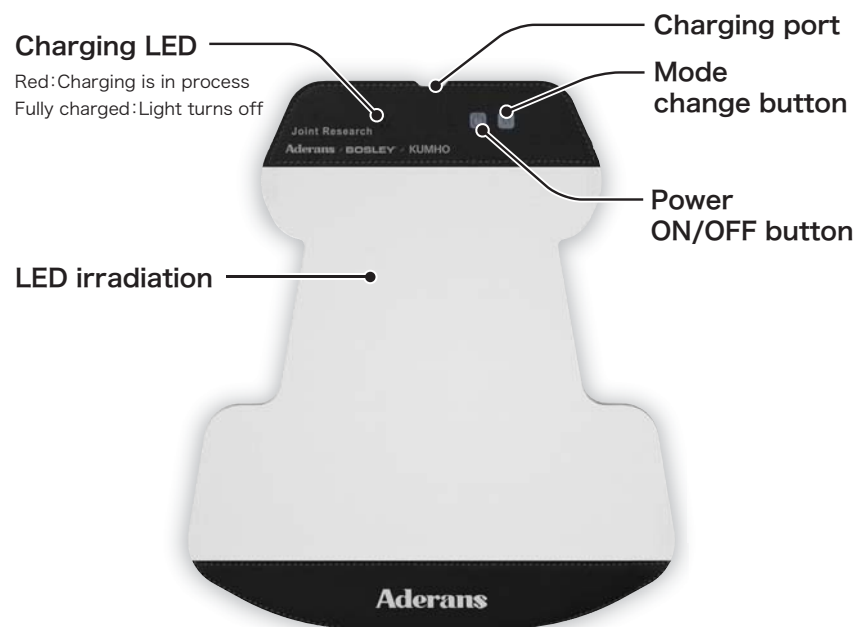Charging LED

Red: Charging is in process
Fully charged: Light turns off

Charging port

Mode change button

Power ON/OFF button

LED irradiation

※This product is a beauty treatment equipment for the scalp and hair, and not for medical purposes.

Product Specification


Name of product	N-LED FLEX
Charging time	Approximately 3 hours
Number of usage after full charging	Approx. 20 minutes x 4 times
Rated input	100-240VAC, 50-60HZ, 0.25A
Rated output	5VDC, 1A
Size/Weight	210.95 x 238 x 3mm (connector 7mm) / approx. 82.52g
Country of origin	Made in Korea
Accessories	Main equipment (1), AC adapter (1), USB cable (1), Instruction Manual & Warranty, Velcro

Mode

Mode	Wavelength range	Lighting order
Mode 1	 630nm	 850nm Red↔Colorless(light red) Color changes every 5 seconds
Mode 2	 415nm  630nm	 850nm Blue, Red, Colorless(light red) Color changes every 5 seconds
Mode 3	 415nm, 630nm	 850nm Pink(red and blue)+Colorless(light red) Color changes every 5 seconds

※Near infrared waves at 850nm are colorless and not visible to the naked eye.


How to use

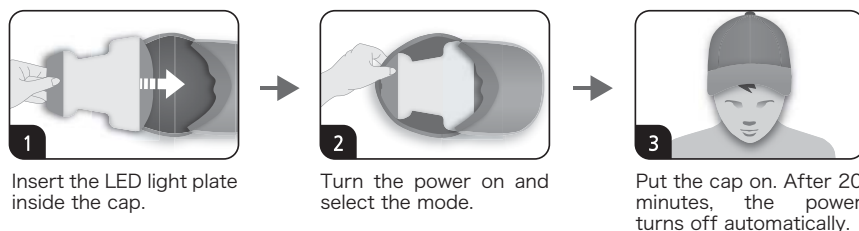 Do not use near the eye area.
Keep out of reach of children.

How to Charge

Red light turns on.
Light turns off when fully charged.

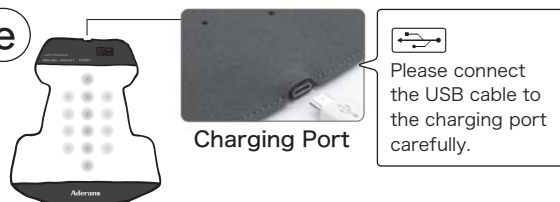Q&A

Q Is it okay to use every day?

A The product is free from side effects on the scalp and hair, even when used every day. However, there can be differences among individuals, thus it is recommended that you constantly check the condition of the scalp and hair to determine the time of use.

A The power turns off automatically after 20 minutes. Take a 10-minute interval before using again. If scalp problems occur, stop using and consult a physician.

Q Is it okay to use while receiving scalp & hair treatment procedures?

A This product is for beauty treatment and not for medical purposes. Follow the doctor's directions according to the condition of the scalp and hair.

Q The charging adapter does not work.

A Check if the product, USB cable and AC adapter are properly connected.

A May not be charged enough. Please try again after charging. Please use it after charging more than 40 minutes.

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